



Install of Media drying H6 and Stand

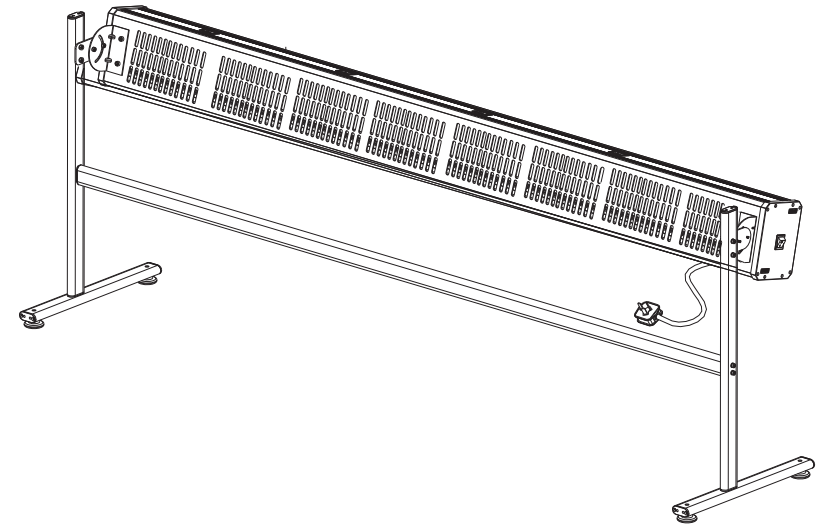

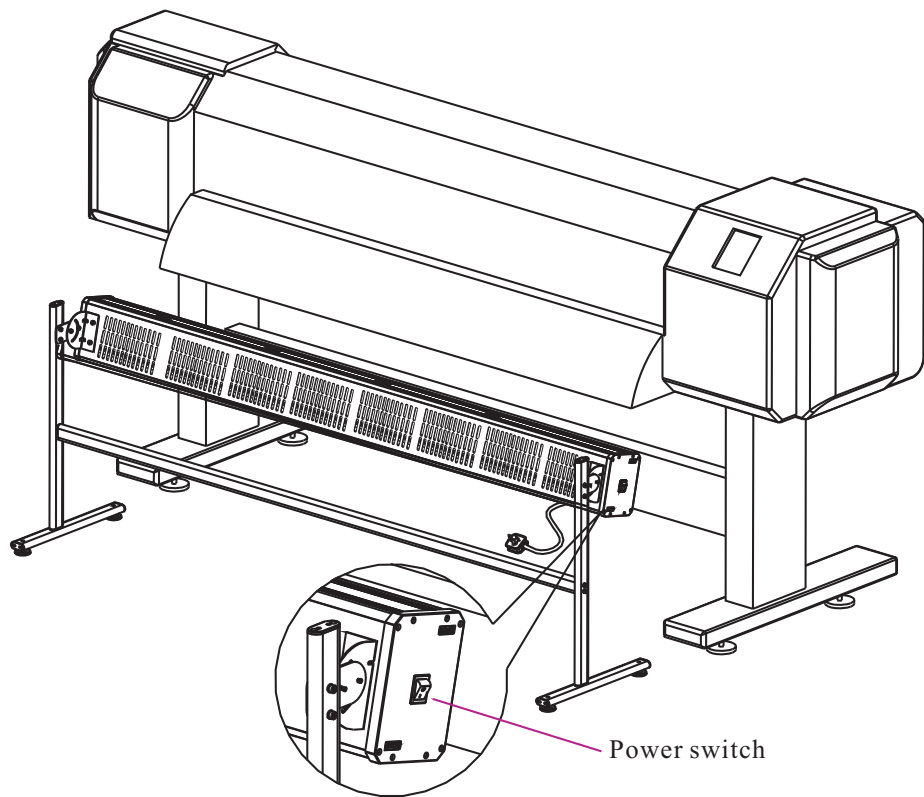

Power switch


Operation Steps:

1. Insert the power plug, Turn on the power switch
2. Turn on the control switch.

Note: Be sure to turn on the power switch, then turn on control switch

Install of Stand


Parts List of Stand

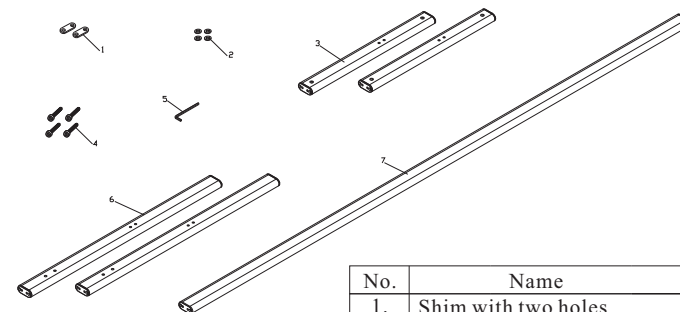

No.	Name	QTY
1.	Shim with two holes	2
2.	Shin \varnothing 5*10	4
3.	Section 1	2
4.	Inner Hexagon Screw M5*35	4
5.	Hexagon Wrench 4mm	1
6.	Section 2	2
7.	Section 3	1


Install drawing of Stand

No.	Name
1.	Inner Hexagon Screw M5*35
2.	shin \varnothing 5*10
3.	shim with two holes

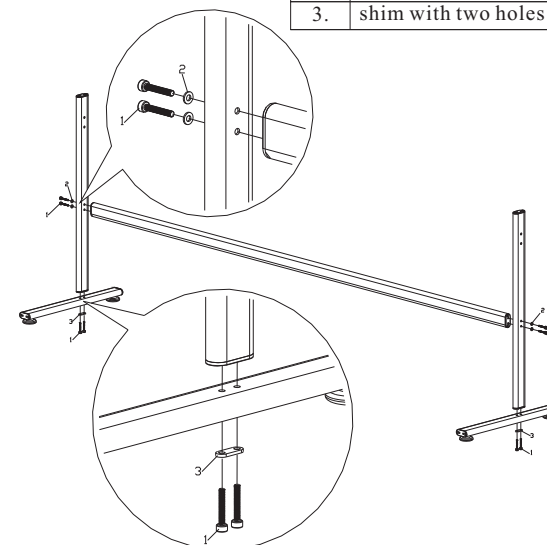

Install of Media drying H6

Parts list of Media drying H6

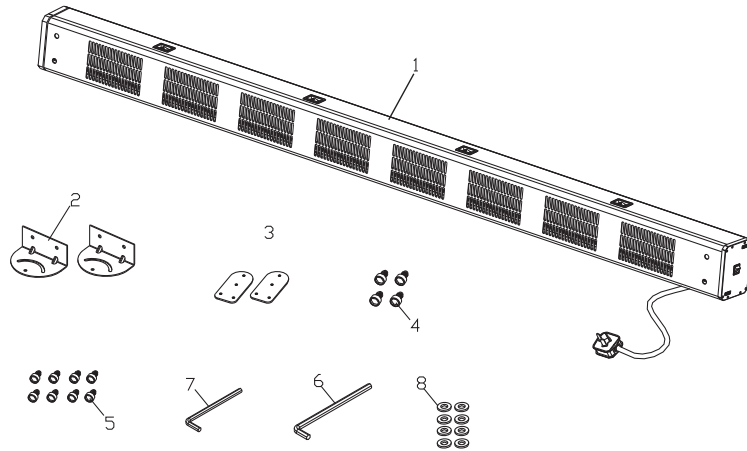


No.	Name	QTY	No.	Name	QTY
1.	Host of Media drying H6	1.	6.	Hexagon wrench 4mm	1.
2.	Slab	2.	7.	Hexagon wrench 5mm	1.
3.	Connection arm	2.	8.	Shin $\varnothing 5*10$	8.
4.	Inner hexagon screw M6*16	4.			
5.	Inner hexagon screw M5*10	8.			

Install of Media drying

Install of Slab and Host of Media drying.
Install requests: at least two persons.


No.	Name
1.	Host of media drying H6
2.	Slab
3.	Inner hexagon screw M6*16

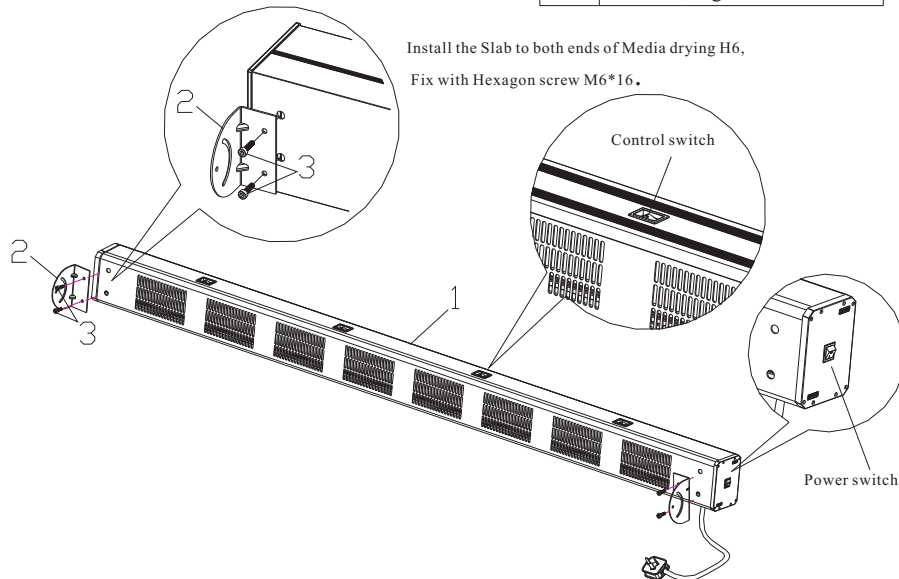

Install of Connecting arm and Slab

No.	Name
1.	Connecting arm
2.	Inner hexagon screw M5*10
3.	Shim $\varnothing 5*10$

Install the Connecting arm to the slab at both ends,
Fix with Hexagon screw M5*10

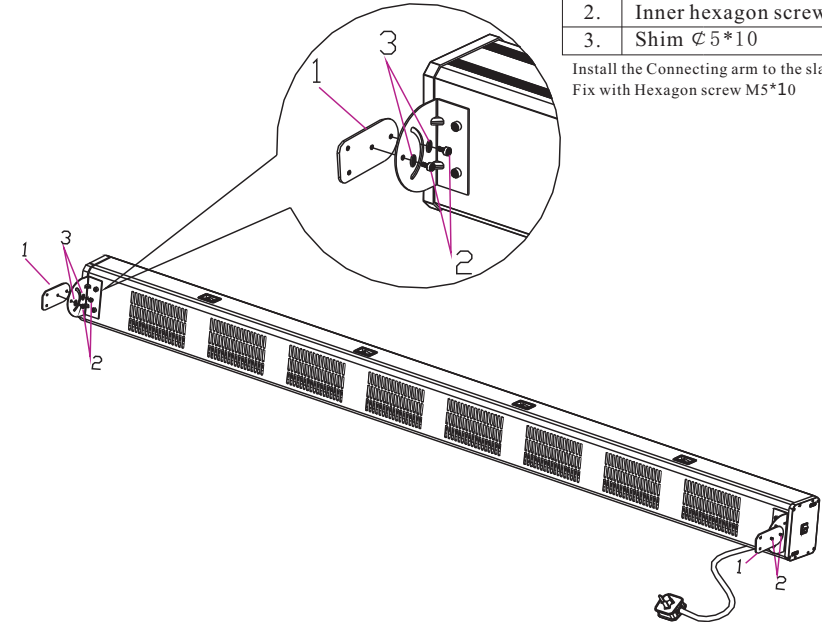

Install of Media drying and Stand

No.	Name
1.	Host of Media drying H6
2.	Stand
3.	Inner hexagon Screw M5*30
4.	Shin $\varnothing 5*10$
5.	Screw Nut M5

